

Heat for infusions
and blood products

 **S-line**



ARTG Number: 167241
Class: 11A

The Barkey *S-line* smart – simple – safe

Advanced technology: One-touch-control

- The Barkey S-line warms blood products and infusions up to the patient.
- The electrically heated and single-sided opening heating device is capable of accepting standard IV-sets.
- The Barkey S-line is simple to operate in all medical departments.
- It is extremely easy to clean due to its smooth surface.
- Its large and flexible application areas at low to medium flow rates result in an optimum alternative to cost-intensive blood heaters with single-use articles.

Avoidance of hypothermia

Warmth is the key to your patients' well-being! With Barkey's inexpensive S-line devices you can apply warmth which eases treatment for all patients and accelerates the healing process. The result: care personnel has less work and treatment costs are reduced.


S-line Made-to-order technology for medicine

Inovative Barkey S-line IV Fluid Warmer - Features & Benefits

- Dry temperature regulation system – no water required
- Integrated temperature sensors – guarantees best warming and highest safety
- Complete encasing of IV sets – optimum warming using standard IV sets
- One touch control – simple to operate in all departments
- Visual & audible alarms – enhanced safety
- Visual display of sleeve temperature – better accuracy & security
- Reusable silicone heating sleeve – latex free care & no consumables
- Preset temperature setting – eliminates the need for manual control
- Smooth Silicone Surface – Easy Cleaning

Exclusive innovative technological solutions at the highest quality standards and specialised benefits in practice - this is the basis for the success of the holistic Barkey heating concept. Barkey has earned its excellent reputation with exclusive engineering know-how in medical technology over the last 25 years. The quality and practical suitability of the devices and systems "made in Germany" has established itself in a wide range of applications.

Type	Barkey S-line
Article number	213.10285
Temperature setting	39°C setpoint temperature ex-works
Heating profile length	L: 1.500 mm
Internal diameters available	4.55 mm 6.55 mm
Suitable for infusion piping diameters	of 4.0 to 7.0 mm
Weight	approx. 750 g
Mains supply	100-240V, 50/60Hz
Connection	24 volt DC
Power consumption	60 Watt
Connecting cable	3 m
Safety class	I protective insulation Degree of protection BF external and internal application to patients with insulated application parts
Humidity protection	IP X2
Certification	CE 0123
Barkey is certified to	ISO 9001 ISO 13485
Made in Germany	

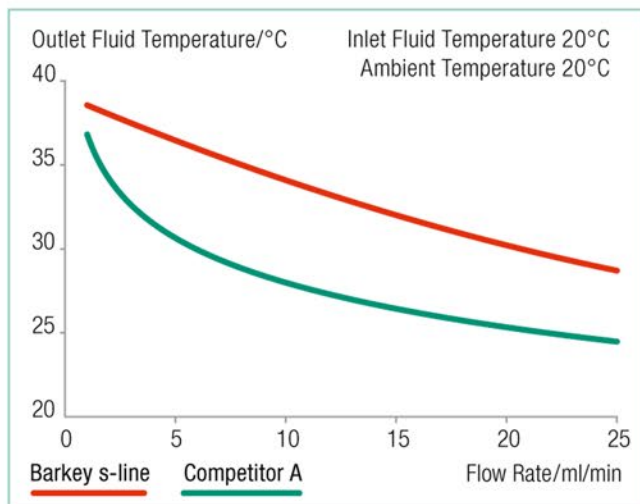

implox Pty.
Ltd.
HEALTHCARE

Units 22-24
60-66 Richmond Road,
Keswick 5035, South Australia
Telephone: 8351 1455
Facsimile: 8293 7377
Email: CustomerService@implox.com
www.implox.com.au